

# F61: NMR–Spectroscopy

T. Gierlich and A. Impertro


May 26th, 2017

# Outline

Introduction and Theoretical Concepts

Part I: Relaxation Times

## Introduction

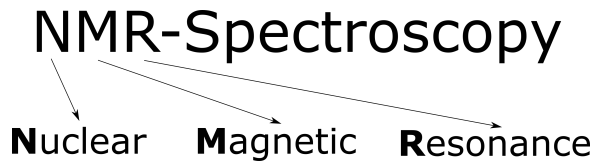

## Introduction


# NMR-Spectroscopy

**Nuclear**

**Magnetic**


**Resonance**

Detection of substances

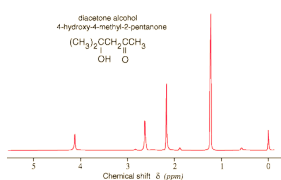


**Figure:** Carl Nave, Hyperphysics, [hyperphysics.phy-astr.gsu.edu](http://hyperphysics.phy-astr.gsu.edu)

Multidimensional imaging

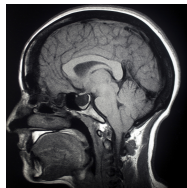

**Figure:** Sierra Vista Diagnostics, [svdrads.com](http://svdrads.com)


# Theoretical Concepts - Working Principle

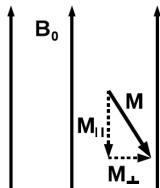

**Figure:** Magnetization


- ▶ Nuclei with spin  $I$  have a magnetic moment  $\mu$
- ▶ Ensemble of many nuclei: Measurable magnetization  $\vec{M}$
- ▶ Minimal energy  $\rightarrow$  Dipole aligned parallel to B-field
- ▶ Ground state  $\rightarrow M_{\perp} = 0$

# Theoretical Concepts - Working Principle

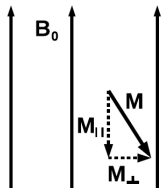

**Figure:** Magnetization

- ▶ Nuclei with spin  $I$  have a magnetic moment  $\mu$
- ▶ Ensemble of many nuclei: Measurable magnetization  $\vec{M}$
- ▶ Minimal energy  $\rightarrow$  Dipole aligned parallel to B-field
- ▶ Ground state  $\rightarrow M_{\perp} = 0$

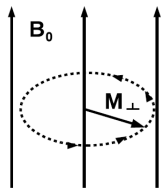

**Figure:** Larmor-Precession of  $M_{\perp}$

- ▶ Excited states have a component  $M_{\perp} \neq 0$
- ▶  $M_{\perp}$  precesses around the field lines with the Larmor frequency

$$\omega_L = \gamma B_0 \quad (1)$$


- ▶  $\omega_L$  can be measured!

# Theoretical Concepts - Working Principle

How can we create an excited state ?

- ▶ An oscillating B-Field  $\vec{B}_1$  rotates the magnetization  $\vec{M}$  by an angle

$$\alpha = \gamma B_1 \Delta t \quad (2)$$

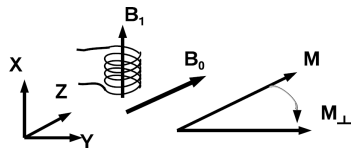


**Figure:** Rotation of  $M$  due to an HF-Pulse

# Theoretical Concepts - Working Principle

## How can we create an excited state ?

- ▶ An oscillating B-Field  $\vec{B}_1$  rotates the magnetization  $\vec{M}$  by an angle

$$\alpha = \gamma B_1 \Delta t \quad (2)$$

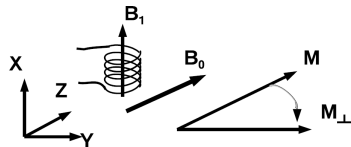

**Figure:** Rotation of  $M$  due to an HF-Pulse

- ▶ By choosing  $\Delta t$ , we can create:
  - ▶ A perpendicular magnetization (90°-Pulse)
  - ▶ An anti-parallel magnetization (180°-Pulse)



# Setup and Measurement Principle

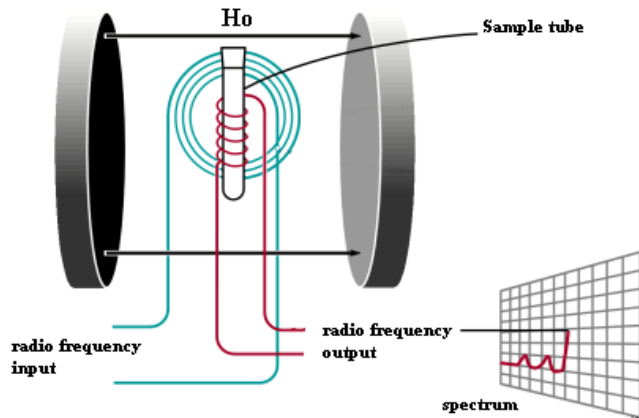


Figure: McGraw Hill Higher Education, mhhe.com

# Theory of Relaxation

**Excited states decay into the Ground State on a characteristic timescale.**  
The decay is of exponential nature and described in the *Bloch equations*:


$$\frac{dM_{\perp}(t)}{dt} = -\frac{M_{\perp}(t)}{T_2} \quad (3)$$

$$\frac{dM_{\parallel}(t)}{dt} = -\frac{M_{\parallel}(t) - M_0}{T_1} \quad (4)$$

- ▶  $T_2$ : Spin-Lattice Relaxation
- ▶  $T_1$ : Spin-Spin Relaxation

# $T_2$ -Measurement: Spin Echo


## Spin-Echo principle

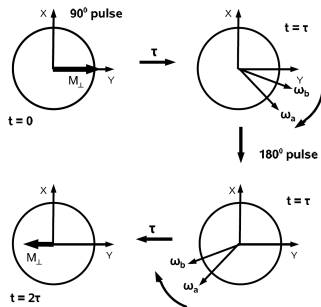

**Figure:** Principle of the spin-echo method


# $T_2$ -Measurement: Spin Echo

## Spin-Echo principle

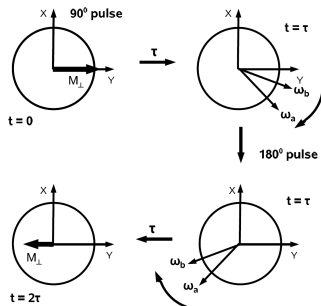

**Figure:** Principle of the spin-echo method

## Pulse sequence

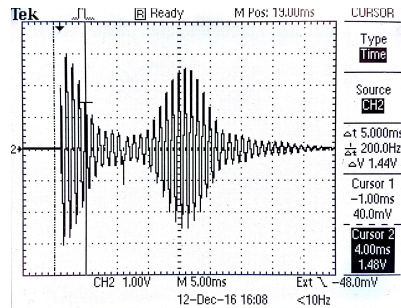

**Figure:** Spin-Echo measurement with  $\tau = 10\text{ms}$

# $T_2$ -Measurement: Spin Echo

## Spin-Echo principle

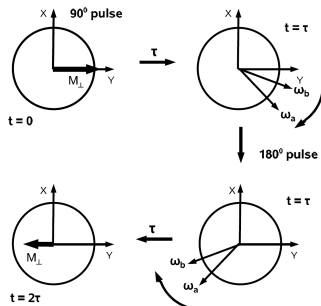


Figure: Principle of the spin-echo method

## Pulse sequence


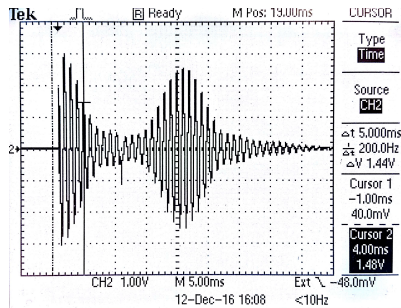



Figure: Spin-Echo measurement with  $\tau = 10\text{ms}$

- Disadvantage: Dephasing for long measurement times!

# $T_2$ -Measurement: Spin Echo

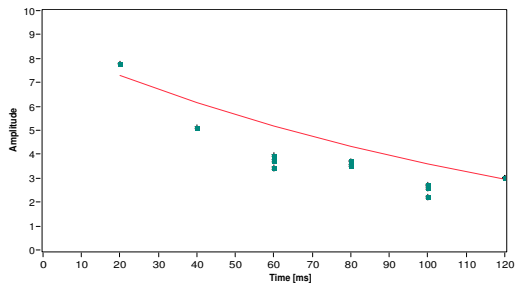

**Figure:** T2-Measurement Sample 1 with fit.

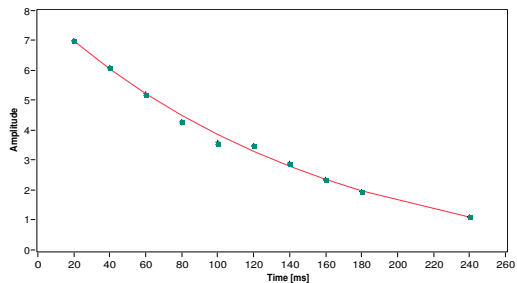

**Figure:** T2-Measurement Sample 3 with fit.

## $T_2$ -Measurement: Carr-Purcell Sequence


**Improve dephasing problem of spin-echo method:**

- ▶ Inject a  $180^\circ$ -Pulse on odd multiples of a time  $\tau$ .
- ▶ The system is phase coherent on even multiples of a time  $\tau$ .


## $T_2$ -Measurement: Carr-Purcell Sequence

Improve dephasing problem of spin-echo method:

- ▶ Inject a  $180^\circ$ -Pulse on odd multiples of a time  $\tau$ .
- ▶ The system is phase coherent on even multiples of a time  $\tau$ .

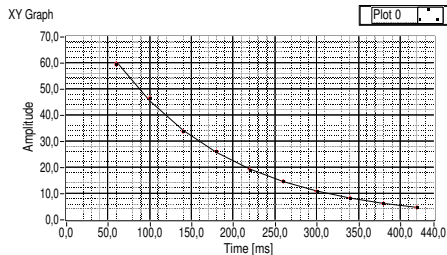

**Figure:** T2-Measurement using Carr-Purcell, Sample 1, with fit.

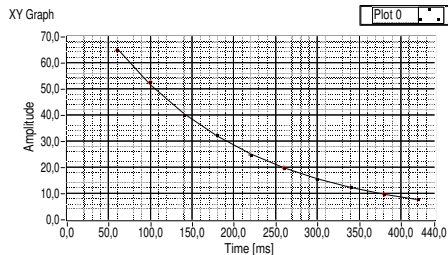

**Figure:** T2-Measurement using Carr-Purcell, Sample 3, with fit.




## $T_1$ -Measurement: Spin Echo


Spin-Echo, but start with a  $180^\circ$ -Pulse (Anti-parallel Magnetization)

# $T_1$ -Measurement: Spin Echo

Spin-Echo, but start with a  $180^\circ$ -Pulse (Anti-parallel Magnetization)

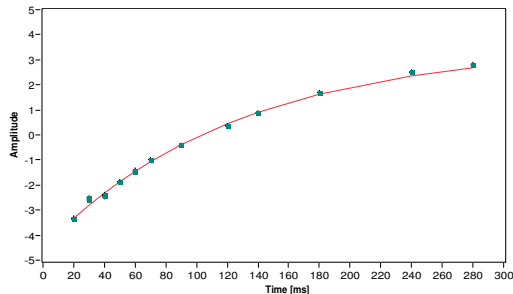

**Figure:** T1-Measurement Sample 1 with fit.

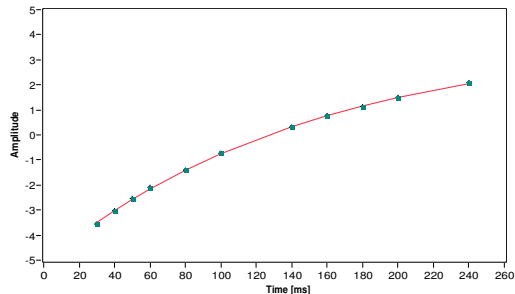

**Figure:** T1-Measurement Sample 3 with fit.

# Relaxation Times: Evaluation

**Table:** Relaxation times – Measured values

Time Method	$T_1$ [ms] Spin-Echo	$T_2$ [ms] Spin-Echo	$T_2$ [ms] Carr-Purcell
Sample 1 (Gd 1:500)	$(125,5 \pm 0,6)$	$(119,5 \pm 0,5)$	$(140,1 \pm 0,4)$
Sample 3 (Gd 1:600)	$(150,5 \pm 1,2)$	$(139,3 \pm 0,8)$	$(166,9 \pm 0,4)$